



LABORATORY SECTOR

SAFETY STORAGE CABINETS



Kalstein®

A different accompaniment, at your service





Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



KALSTEIN® Safety Storage Cabinets
2-4 Gallon Explosion Proof Cabinet
Model YR05464 // YR05465



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® Safety Storage Cabinets 2-4 Gallon Explosion Proof Cabinet Model YR05464 // YR05465





Product Description:

- ✓ Environmental protection spray, strong insulation
- ✓ Three-point linkage, double lock
- ✓ Strong load-bearing, adjustable laminate
- ✓ Electrostatic grounding, safety upgrade
- ✓ Three colors, each with its own strengths

Features:

- ✓ Double lock safe and reliable: Two locks, two managers. Double security, pad lock and cabinet lock, standardized management, improve safety factor.
- ✓ Three point linkage solid lock post: Three point linkage cabinet door lock tongue locking structure, upper, middle and lower three-point linkage solid core lock columns to increase the security of the cabinet.
- ✓ Adjustable fire vent: The fire vents on both sides of the cabinet body are equipped with adjustable valves, which can rotate to control the air volume and reduce the risk of harmful gas accumulation in the cabinet.
- ✓ Liquid leak proof galvanized shelves: 4-110 gallons, regardless of size, the shelves are designed with leaked liquid guide groove, which is not afraid of leakage. The laminate can be adjusted up and down according to demand.
- ✓ Anti-static: each model is equipped with OSHA standard electrostatic prevention wire clip. Ensure safe use and management.

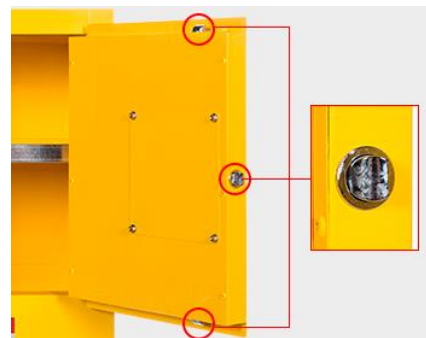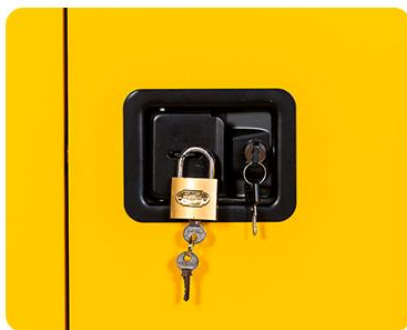

#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® Safety Storage Cabinets
2-4 Gallon Explosion Proof Cabinet
Model YR05464 // YR05465



Technical Description:

| Model | YR05464 | YR05464-1 | YR05464-2 | YR05465 | YR05465-1 | YR05465-2 |
|-------------------------|------------------------------|-----------------------------|----------------------------|------------------------------|-----------------------------|----------------------------|
| Dimensions | 330 x 330 x 330 | | | 560 x 430 x 430 | | |
| (H x W x D) | | | | | | |
| Weight (Kg) | | | | 21 | | |
| Volume | 02 GAL / 7 L | | | 04 GAL / 15 L | | |
| Door Type | Single door / Manual | | | Single door / Manual | | |
| Description | Yellow for Flammable Liquids | Red for Combustible liquids | Blue for Corrosive liquids | Yellow for Flammable Liquids | Red for Combustible liquids | Blue for Corrosive liquids |
| Adjustable shelf | 1 Layer | | | | | |



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
 2 Rue Jean Lantier • 75001 Paris •
 +33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S





 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



**Fundación
Jacinto Convit**

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



Humatem

Solidarité et coopération
biomédicale


<https://www.humatem.org/>



**MANIAPURE
FOUNDATION USA**

<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)

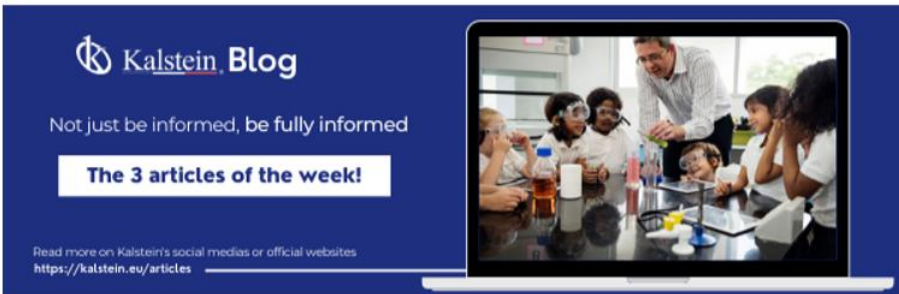

Kalstein Blog
Not just be informed, be fully informed

The 3 articles of the week!

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



1
Read this ARTICLE NOW

How does the climate affect your microscope?

Kalstein Blog

Article of the week:
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW

What is the importance of the laboratory refrigerator?

Kalstein Blog

Article of the week:
What is the importance of the laboratory refrigerator?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...

23 views • 3 years ago



Kalstein's Microscopes

445 views • 1 week ago



PARIS - FRANCE
HEADQUARTER

2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89
Cel: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

