




LABORATORY LINES

BACTI-CINERATOR STERILIZER



Kalstein®

A different accompaniment, at your service

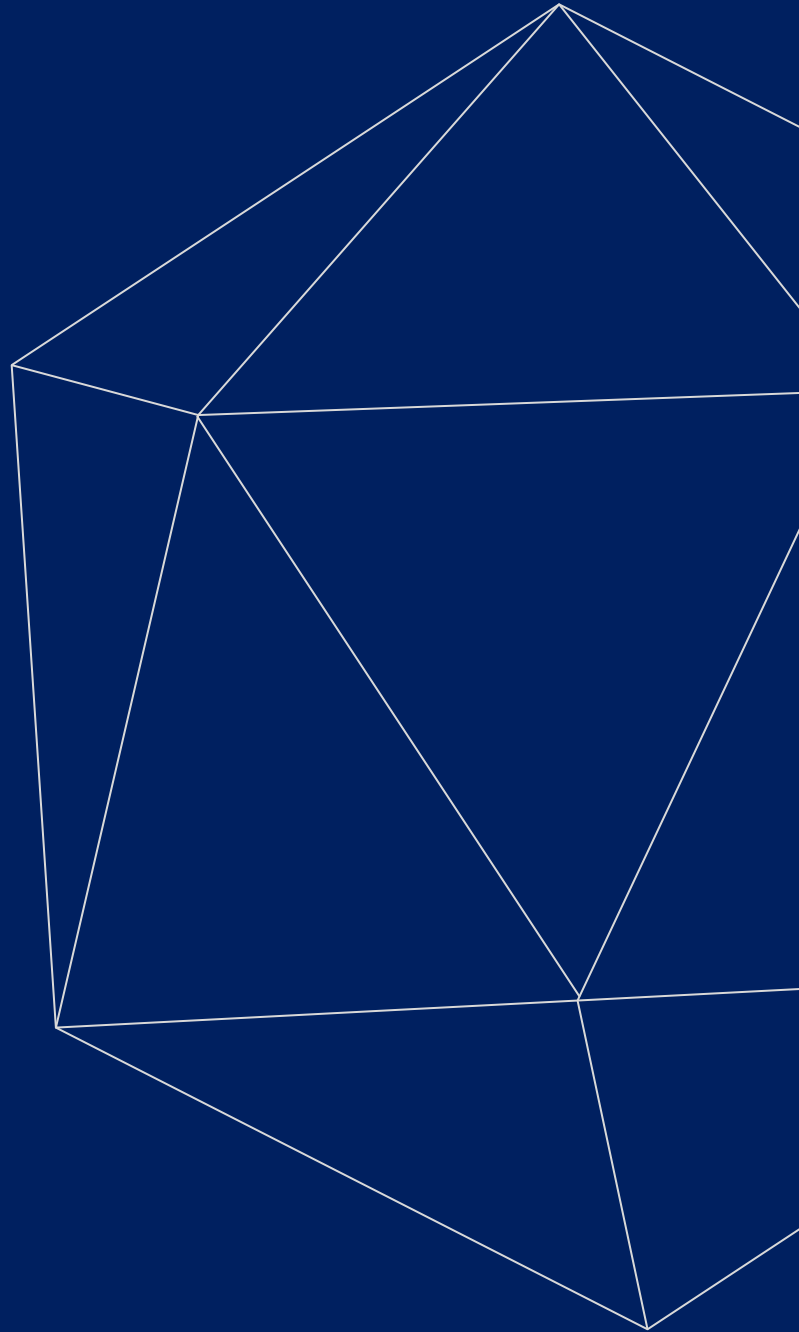



Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



KALSTEIN® Bacti-cinerator Sterilizer

YR05324 - YR05325



Kalstein's bacti-cinerator sterilizer sterilizes micro-organisms utilizing infrared heat produced by a ceramic core element. The ceramic element contains no asbestos and ensures maximum sterility without splatter across the work surface area. Complete sterilization occurs within 5-7 seconds at optimum sterilizing temperature of 1500°F (815.6°C).

The small footprint makes the bacti-cinerator sterilizer an ideal instrument for anaerobic and aerobic chambers alike. Several user friendly accessories are also available for the bacti-cinerator sterilizer.

The YR05324 and YR05325 are our company new products for the bacti-cinerator sterilizer. The angle of the heater will be changed, and it's very convenient for sterilizing the large flask, bottle and etc.

Features

- Sterilizes needles, loops, and culture tube mouths in five to seven seconds at optimum sterilizing temperature of 815.6°C (1500°F).
- Ceramic funnel tube enclosed in stainless-steel perforated guard and stainless-steel support stand.
- For AC 220V/AC120V 50/60Hz.



YR05324



YR05325

Technical Specifications

Model	YR05324	YR05325
Center high-temperature	825°C ± 50°C	825°C ± 50°C
Max-diameter of sterilizer goods	Φ14mm	Φ35mm
Length of sterilizer	150mm	100mm
Dimensions (DxWxH)mm	159x100x260	153x100x280
Power	170W	320W
Weight	1.5kg	1.6kg
Adjustable angle range	55°	55°



#Let'sgivemore ❤️

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S





A different accompaniment, at your service

WITH THE ACQUISITION OF A KALSTEIN EQUIPMENT

YOU MAKE A CONTRIBUTION TO:

#LETSGIVEMORE



Fundación Jacinto Convit

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



Humatem

Solidarité et coopération biomédicale

<https://www.humatem.org/>




MANIAPURE FOUNDATION USA

<https://www.maniapurefoundation.org/>

Thank you!

Articles Kalstein
(Social Networks and Official Websites)

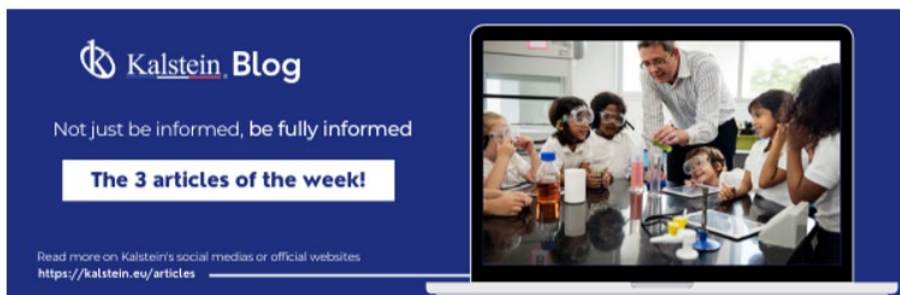

Kalstein Blog
Not just be informed, be fully informed
The 3 articles of the week!
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



1
Read this ARTICLE NOW
How does the climate affect your microscope?
A microscope is an instrument that allows observing and measuring objects at a very small scale. This is achieved through a complex system of lenses and mirrors, which form...

Kalstein Blog
Article of the week:
How does the climate affect your microscope?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW
What is the importance of the laboratory refrigerator?
A refrigerator is a device that is one of the most common pieces of equipment used in laboratories. It is a controlled environment (refrigerated space), various fluids and substances...

Kalstein Blog
Article of the week:
What is the importance of the laboratory refrigerator?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Alliances
23 views • 3 years ago



Kalstein's Microscopes
445 views • 1 week ago



PARIS - FRANCE
HEADQUARTER
2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89
Cel: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

