




LABORATORY LINES

ICE MAKER



Kalstein®

A different accompaniment, at your service

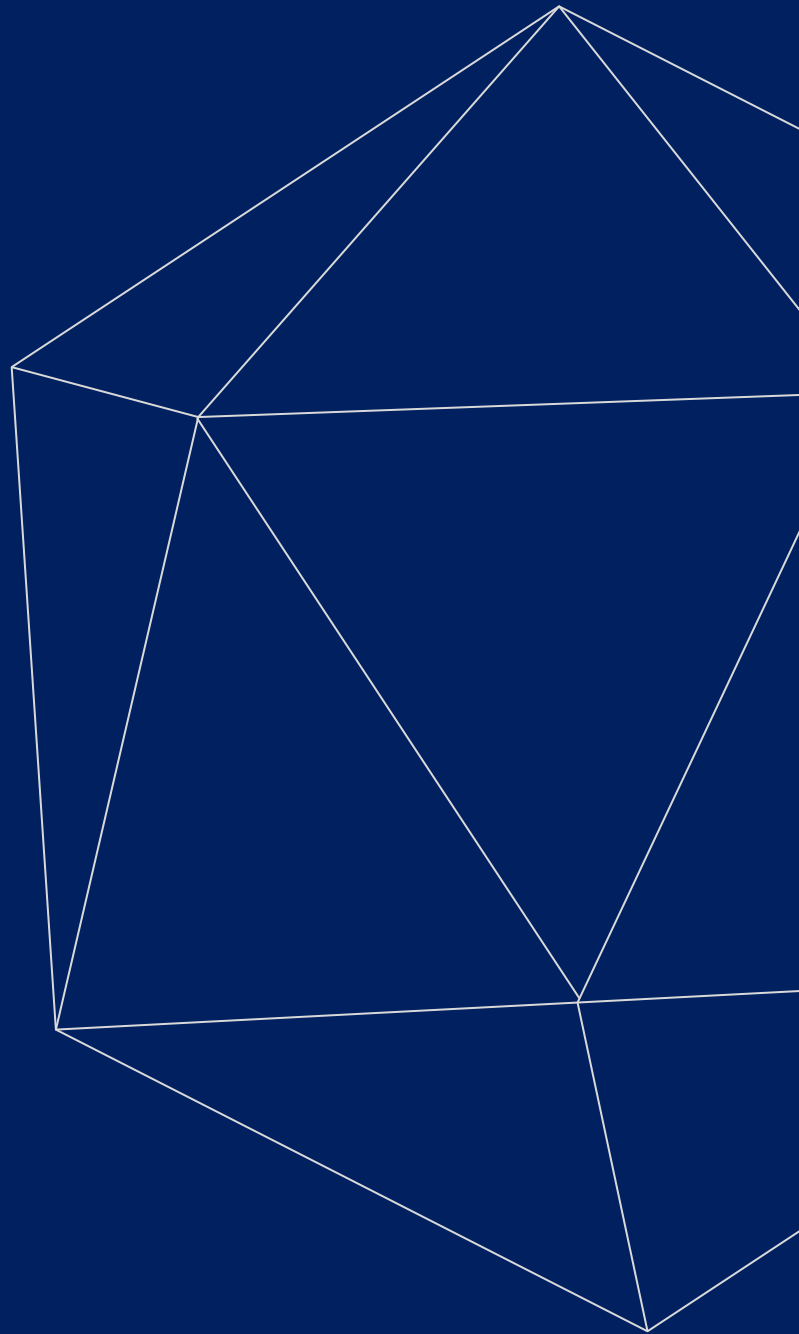



Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore ♥

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.





Flake Ice Maker Machine

Serie YR05134/ YR05135/YR05136

Features:

- ✓ High grade stainless steel is used which are anti-corrosive and durable. The compact structure is simple and space-saving.
- ✓ The cabinet are filled with fluoride-free foam thermal cover featuring good insulation effect. The inner liner is food grade plastic. They are both energy-saving and environment-friendly.
- ✓ World famous brand CFC free R134a compressors are used. Major parts and components have passed relevant safety certifications. The product has stable quality and long service life.
- ✓ The ice-making process is fully computer controlled, using imported computer chips to make control reliable and operation easily.
- ✓ The advanced technologies and imported accessories makes noise low, operation smooth and reliable. The ventilation holes and fans are deployed to ensure that motor of reducer can work reliably with high ambient temperature.
- ✓ The patented technology of line of chip ice-making chamber separated from the evaporator which provides the high cooling efficiency and large ice-making capacity.
- ✓ The screw extrusion hob-style ice-making type makes a compact structure, and automatic separation of ice from water. The optimal design of skate blade makes the shape of ice small and practical.
- ✓ The unique float-type water tank water inflow system ensures that there is no residual water. Ice clearance process not needed, and no water wasted, saving water and energy.
- ✓ Fully automatic running with ice full display, water shortage display, in operation display and malfunction warning display, etc. for shutdown protection. They will stop automatically in case of ice full or water shortage and start automatically when ice is not full or water supply is resumed.
- ✓ Flake ice cools more quickly than other ice forms, which can infiltrate into narrow space, featuring fast cooling speed, good ice bathing effect.
- ✓ Long-term performance tests and safety electrical tests will be performed before delivery, ensuring excellent performance of the product.



#Letsgivingmore 

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



Flake Ice Maker Machine

Serie YR05134/ YR05135/YR05136



Parameters per Model

Model:	YR05134	YR05135	YR05136
Ice Making Capacity:	60Kg/24h	70Kg/24h	85Kg/24h
Ice Storage :	25Kg	25Kg	25Kg
Voltage:	220-240V/50Hz	220-240V/50Hz	220-240V/50Hz
Power:	420W	460W	480W
Cooling type:	Air Cooling	Air Cooling	Air Cooling
Refrigerant:	R134a	R134a	R134a
Water Supply:	Tap Water	Tap Water	Tap Water
Net Weight: Kg	42	42	42
Gross Weight: Kg(Carton)	45	45	45
Product Dimension: (mm)	400*510*840	400*510*840	400*510*840 mm
Carton Packaging: (mm)	470*580*890	470*580*890	470*580*890 mm



#Letsgivemore ❤️

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
 2 Rue Jean Lantier • 75001 Paris •
 +33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S.



Flake Ice Maker Machine

Serie YR05134/
YR05135/YR05136

01

Thickened Insulation Layer


Strong insulation bin, ice can be stored for 5-8 hours after ice full shutdown.



02

SUS304

High-quality 304 stainless steel, beautiful and durable.

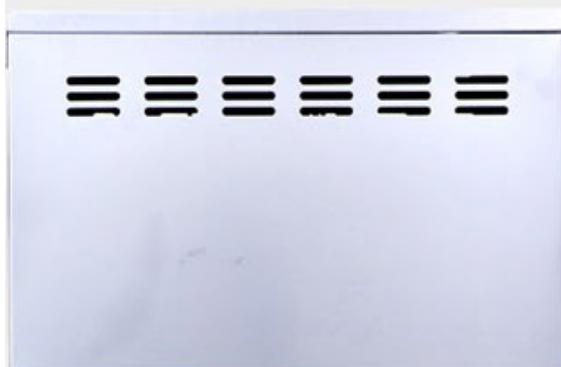

03

Water Filter

Set the filter at the inlet, to prevent large particles into the machine.



04

All-round Cooling

Cooling holes are installed around the machine, to avoid high temperature operation of the compressor.



#Letsgivemore 

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S





Flake Ice Maker Machine

Serie YR05134/ YR05135/YR05136



#Letsgivingmore 

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S





 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



**Fundación
Jacinto Convit**

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



Humatem

Solidarité et coopération
biomédicale


<https://www.humatem.org/>



**MANIAPURE
FOUNDATION USA**

<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)

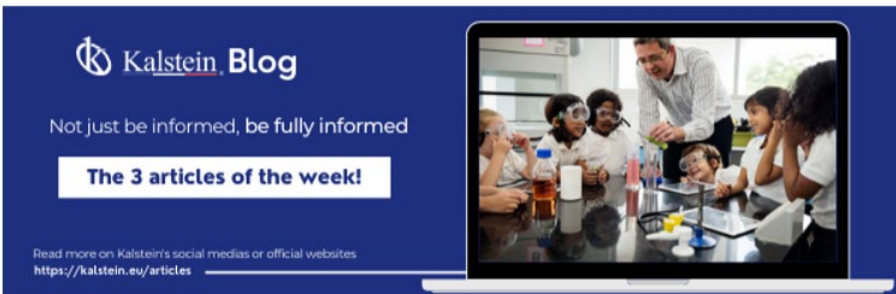

Kalstein Blog
Not just be informed, be fully informed

The 3 articles of the week!

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

1
Read this ARTICLE NOW



How does the climate affect your microscope?

microscope is an instrument that allows observing and measuring objects at a very small scale. This is achieved through a complex system of lenses and mirrors, which form

Kalstein Blog

Article of the week:
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

2
Read this ARTICLE NOW



What is the importance of the laboratory refrigerator?

refrigerator is a device that is one of the most common pieces of equipment used in laboratories. It is a controlled environment (refrigerated space), various fluids and substances

Kalstein Blog

Article of the week:
What is the importance of the laboratory refrigerator?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Alliances
23 views • 3 years ago



Kalstein's Microscopes
445 views • 1 week ago



PARIS - FRANCE
HEADQUARTER
2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89
Cel: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

