



LABORATORY SECTOR

ICE MAKER



Kalstein®

A different accompaniment, at your service





Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



KALSTEIN® Ice Maker
High Quality Shaved Ice Machine For Laboratory
Model YR05140 // YR05141



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® Ice Maker

High Quality Shaved Ice Machine For Laboratory

Model YR05140 // YR05141



Product Description:

- ✓ High-quality stainless steel is used anti-corrosion and durable. The compact structure is simple and saves space.
- ✓ The cabinet is filled with fluorine-free thermal foam with good insulation effect. The inner lining is food grade plastic. They save energy and are respectful with the environment.
- ✓ World famous brand CFC-free R134a compressors are used. The main parts and components have passed the relevant safety certifications. The product has stable quality and long service life.
- ✓ The ice making process is fully computer controlled, using imported computer chips to make the control reliable and easy to operate.
- ✓ Advanced technologies and imported accessories make low noise, smooth and reliable operation. The vents and fans are deployed to ensure that the geared motor operates reliably.
- ✓ The patented technology of the Chip Ice Chamber line separated from the evaporator that provides high cooling efficiency and large ice making capacity.
- ✓ Screw extrusion ice making type makes compact structure and automatic separation of ice from water. The optimal design of the skate blade makes the ice shape small and practical.
- ✓ Unique float type water inlet system ensures no residual water. no need for an ice cleaning process or wasted water, saving water and energy.
- ✓ Fully automatic operation with full ice screen, water shortage screen, operation screen and malfunction warning screen, etc. For power-off protection. They will stop automatically in case of ice or water shortage and start automatically when the ice is not full or the water supply is restored.
- ✓ Flake ice cools faster than other forms of ice, which can infiltrate a narrow space, with fast cooling speed, good ice bath effect.
- ✓ Long-term performance testing and electrical safety testing will be conducted before shipment, ensuring excellent product performance



#Letsgivemore 

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® Ice Maker High Quality Shaved Ice Machine For Laboratory Model YR05140 // YR05141



Technical Description:

Model		YR05140	YR05141
Ice Making Capacity		200 Kg / 24 h	250 Kg / 24 h
Ice Storage		60 Kg	
Voltage		220V / 50Hz or 110V / 60HZ (you need to pay extra)	
Power		1040 W	1100 W
Cooling Type		Air cooling	
Refrigerant		R134a	
Water Supply		Tap water	
Net Weight (Kg)		80	
Gross Weight (Kg)		85	
Product Dimension		680 x 611 x 1 100 mm	
Packaging Dimension		880 x 750 x 1 250 mm	
Filling quantity	20'GP	48	
	40'GP	102	
	40'HQ	102	



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S





 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**


YOU MAKE A CONTRIBUTION TO:



<https://onetreepanted.org>



Articles Kalstein
(Social Networks and Official Websites)

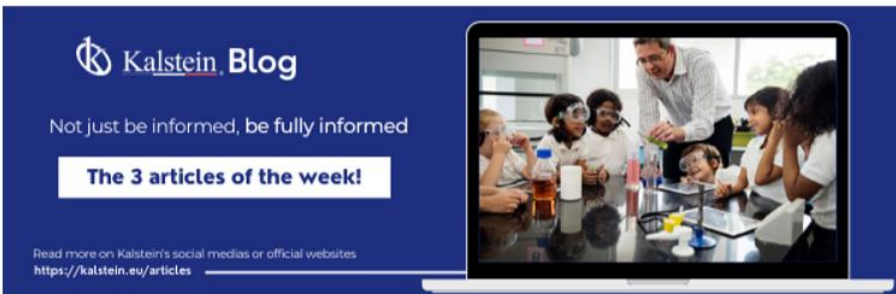

Kalstein Blog
Not just be informed, be fully informed

The 3 articles of the week!

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

1
Read this ARTICLE NOW



How does the climate affect your microscope?

Kalstein Blog

Article of the week:
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

2
Read this ARTICLE NOW



What is the importance of the laboratory refrigerator?

Kalstein Blog

Article of the week:
What is the importance of the laboratory refrigerator?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...
23 views • 3 years ago



Kalstein's Microscopes
445 views • 1 week ago



PARIS - FRANCE
HEADQUARTER
2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89
Cel: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

