




MEDICAL LINES

INFANT WARMER



Kalstein®

A different accompaniment, at your service

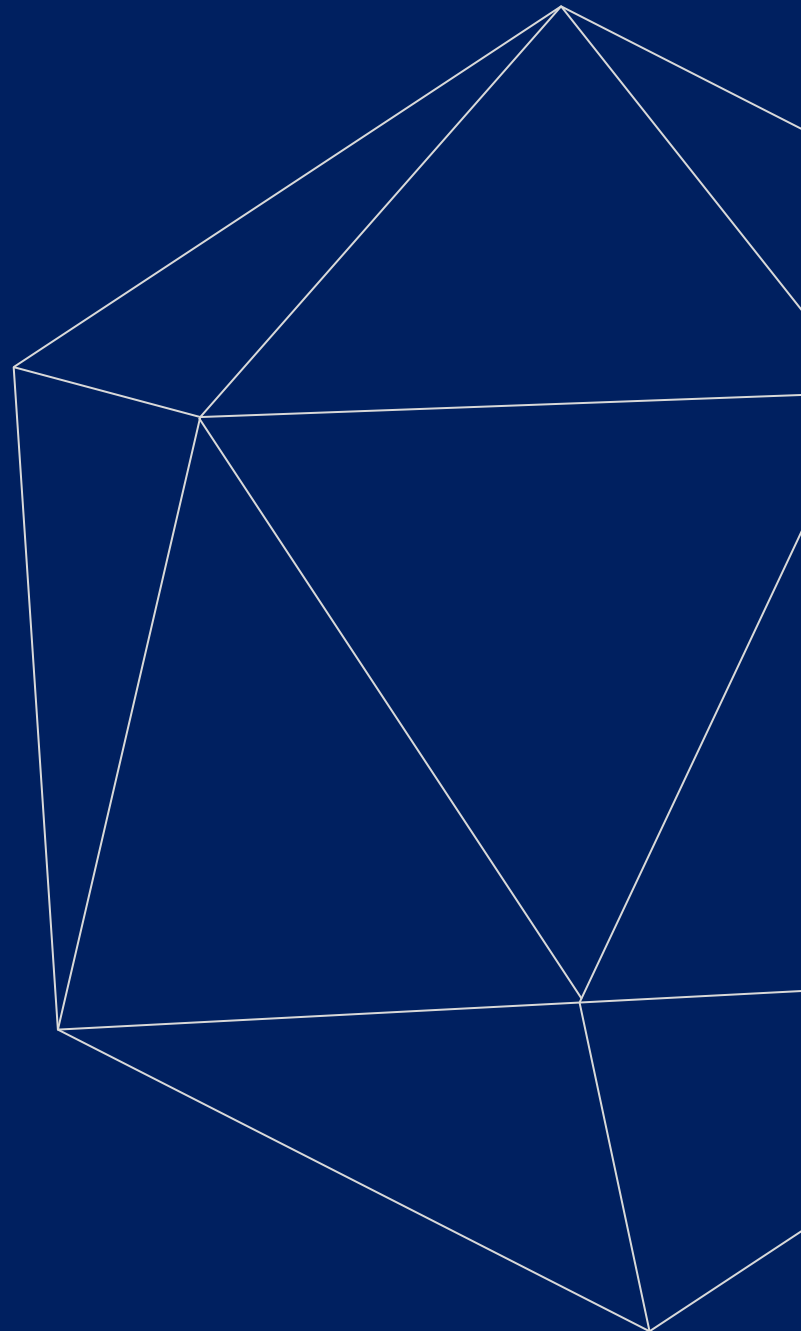



Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



KALSTEIN® Infant Radiant Warmer

Infant Radiant Warmer

Special Features



YR02191



YR02189



YR02190

Parameter system

	YR02189	YR02190.....	YR02191	YR02192
Temp control mode	*Pre-warming, Manual control and Skin temperature*			
Rotatable radiant warmer	●	●	●	●
Tilttable bed	●	●	●	●
Failure alarm	●	●	●	●
APGAR timer	○	●	●	●
X-ray cassette	○	●	●	●
Bed dimension	standard	standard	large	large
Wings of phototherapy units	○	○	○	○
Resuscitation bag	○	○	○	○
Adjustable working height	N/A	N/A	○	●
Low- pressure suction unit	○	○	○	○
Oxygen therapy unit	○	○	○	○
Oxygen hood	○	○	○	○
Infant head restraint unit	○	○	○	○

"○" means optional "●" means standard configuration



#Letsgivingmore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.




All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



KALSTEIN® Infant Radiant Warmer

Infant Radiant Warmer

Special Features

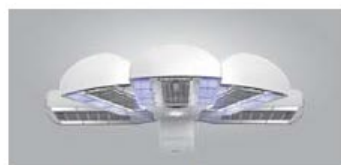

Wing-type phototherapy units for all model (optional)



Compact swiveled shelf for taking monitors for YR02191 only



X type base with shelf for storage design for YR02191 only



Vertical height adjustment with ergonomic design for YR02192 only



Basic base with shelf for storage design for YR02189/ YR02190



The warmer module can be swiveled from $\pm 90^\circ$ for taking X-rays for YR02191/YR02192 / YR02190



The warmer module can be swiveled from $\pm 90^\circ$ for taking X-rays for YR02189



Baby bed can be adjusted up to $\pm 15^\circ$ with gas spring design for YR02191/YR02192



X-ray tray beneath mattress for your convenience of tacking X-ray for YR02192/ YR02191/ YR02190



Bigger size warming bed is enclosed by four protective panels for YR02191/ YR02192

Features

- ✓ Temp. control modes: Pre-warming, Manual control and Skin temperature
- ✓ Both carborundum & far-IR electric or crystal quartz heater available
- ✓ Adjustable operating height with lift system
- ✓ APGAR score & timer system for evaluating babies' Activity, Pulse, Grimace, Appearance and Respiration situation
- ✓ X-ray cassette tray under bed for easy imaging
- ✓ Dual-channel microprocessor controlled
- ✓ Temperature deviation observation and amendement on front panel
- ✓ Anti-separation function of skin temperature sensor for over-warming prevention
- ✓ RS-232 interface for networking and data communication



#Let'sgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



KALSTEIN® Infant Radiant Warmer

Infant Radiant Warmer

Model	YR02193	YR02192	YR02191	YR02190	YR02189
General					
Control system	Servo & Manual Control by Microprocessor Control System				
Display	LED digital display				
Mattress size	57×76cm	57×76cm	57×76cm	57×76cm	54×66cm
Baby bed tilt	±15°continuous tilt	±15°continuous tilt	±15°continuous tilt	±15°continuous tilt	±10°continuous tilt
Power supply	110 ~ 120V / 220 ~ 240V, 50/60 Hz				
Heater source					
Far-infrared heater	silicon carbide	silicon carbide	silicon carbide	silicon carbide	anti-explosion crystal quartz
Consumption	≤800VA				
Heater output indicator	0-100%, 10% increments.				
Horizontal angle adjustment of heater lamp	±90°				
Temperature control mode					
Control mode	Preheating, manual, skin-temp mode				
Temperature display range	0.0 °C~99.0 °C				
Temperature resolution	0.1 °C				
Skin-temp control mode					
Skin-temp control range	25.0 °C~37.0 °C				
Temperature display range	25.0 °C~42 °C				
Temperature resolution	0.1 °C				
Wings of infant phototherapy unit					
Light source	LED light	LED light	LED light	optional	optional
Power output	≤200W	≤200W	≤200W	optional	optional
Lifetime	20000hrs	20000hrs	20000hrs	optional	optional
Built-in timer range	0~9999hrs 59 min	0~9999hrs 59 min	0~9999hrs 59 min	optional	optional
Wavelength	460nm	460nm	460nm	optional	optional
Electrical height adjustment	yes	yes	optional	no	no
APGAR timer	1min / 5min / 10min				
Alarm system					
Sensor failure	Over-high temp / low temp				
Temperature deviation	Yes				
Power failure	Yes				
Heater failure	Yes				
System failure	Yes				
Configurations					
Storage drawer	optional				
Baby weight scales	optional				
Radiolucent mattress	optional				
Down side LED phototherapy of baby bed	standard	standard	standard	optional	optional
Resuscitation bag	optional				
Low-pressure suction unit	standard	optional	optional	optional	optional
Oxygen therapy unit	optional				
Oxygen hood	optional				
Infant head fixing unit	optional				



#Let'sgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S





 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



**Fundación
Jacinto Convit**

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



Humatem

Solidarité et coopération
biomédicale


<https://www.humatem.org/>



**MANIAPURE
FOUNDATION USA**

<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)

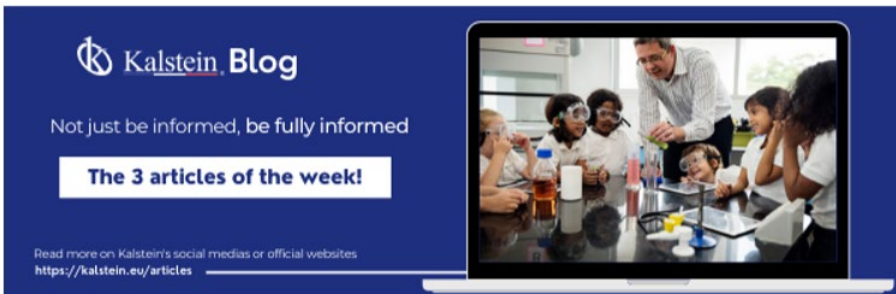

Kalstein Blog
Not just be informed, be fully informed
The 3 articles of the week!
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



1
Read this ARTICLE NOW
How does the climate affect your microscope?
A microscope is an instrument that allows observing and measuring objects that are too small to be seen with the naked eye. This is achieved through a complex system of lenses, which form
Kalstein Blog
Article of the week:
How does the climate affect your microscope?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW
What is the importance of the laboratory refrigerator?
A refrigerator in a laboratory is one of the most important equipment for scientists to maintain, in a controlled environment (refrigerated space), various fluids and substances.
Kalstein Blog
Article of the week:
What is the importance of the laboratory refrigerator?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Alliances

23 views • 3 years ago



Kalstein's Microscopes

445 views • 1 week ago



MONTPELLIER - FRANCE
HEADQUARTER
Optimum Business Center,
450 Rue Baden Powell,
34000 Montpellier - France.

Fax: +33 (0) 4 67 15 88 49
Tel: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

Official Websites for Latin America.

