



LABORATORY SECTOR

SEALING MACHINE



Kalstein®

A different accompaniment, at your service





Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



KALSTEIN® Sealing Machine
Sealing Machine
Model YR05987



Product Description:

Sterilization in an autoclave is the safest way for prevention of infections. After sterilization, however, care must be taken that instruments are not contaminated again. The Sealing Machine YR05987 can effectively seal sterilized instruments so that they are stored while being safely protected against contamination. Without great physical effort, these two models produce optimal sealing seams 10 mm wide, with precise sealing pressure. The German standard DIN 58953.

Specifies sealing seams at least 8 mm wide. The working table is integrated in the housing of the device. After the operator switches on the power, it takes approx. 2 minutes until the cold unit heats up to operating temperature. When the signal lamp goes off, this means that the YR05987 is in the ready mode and can go to work. When the operator turns the crank, this presses the pressure bar against the heated aluminum rail. After only 2 seconds, the sterilization wrapper has been properly sealed in accordance with relevant standards, with a sealing seam of 10 mm width.



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
 2 Rue Jean Lantier • 75001 Paris •
 +33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® Sealing Machine

Sealing Machine

Model YR05987



Technical Data:

Model	YR05987
Voltage	220V/50Hz
Power	500 W
Sealing Width	250mm
Press Width	10mm
Outside size	480×410×290mm
Net Weight	5.5 kg

Features:

- ✓ Maximum Sealing-Sealing temperature adjusts automatically when voltage changes.
- ✓ Non-corrosive body-the device is made of stainless for easy cleaning.
- ✓ Sharp cutting blades-two cutting blades are installed to cut evenly.
- ✓ It is made of sterilization paper and polypropylene, coated with thermal chemical indicator.
- ✓ It is used to bag instruments for sterilization. Instruments will remain sterile for a long time.
- ✓ There are different widths of 5 cm,7 cm and 10 cm to best suit instrument sizes.



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S





 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



**Fundación
Jacinto Convit**

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



Humatem

Solidarité et coopération
biomédicale


<https://www.humatem.org/>



**MANIAPURE
FOUNDATION USA**

<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)

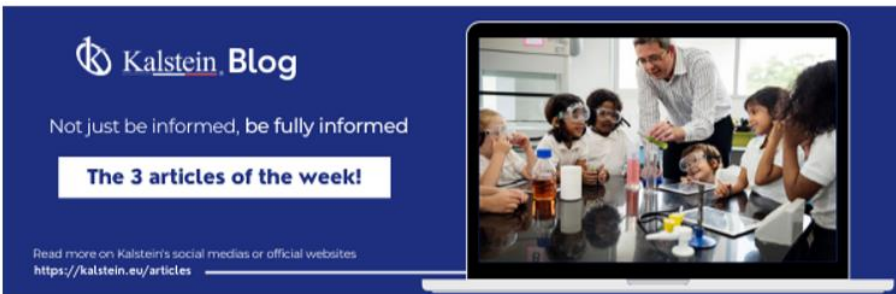

Kalstein Blog
Not just be informed, be fully informed

The 3 articles of the week!

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

1
Read this ARTICLE NOW



How does the climate affect your microscope?

Kalstein Blog

Article of the week:
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

2
Read this ARTICLE NOW



What is the importance of the laboratory refrigerator?

Kalstein Blog


Article of the week:
What is the importance of the laboratory refrigerator?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...
23 views • 3 years ago

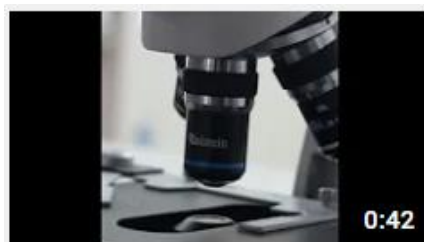

Kalstein's Microscopes
445 views • 1 week ago



PARIS - FRANCE
HEADQUARTER
2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89
Cel: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

