

LABORATORY SECTOR

ULTRASONIC CLEANER







All rights reserved ® KALSTEIN France S. A. S., 2 Rue Jean Lantier • 75001 Paris • +33 1 78 95 87 89 / +33 6 80 76 07 10 • https://kalstein.eu KALSTEIN FRANCE, S.A.S





OUR SERVICES Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings In any part of the world, receive your

induction or training from our specialized team of engineers



Quick Response Our work team is always available to response all your consults or questions, in order to support you in any situation.





#Letsgivemore 💙

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or.air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.





KALSTEIN® Ultrasonic Cleaner Ultrasonic Washing Machine Portable Model YR05207





#Letsgivemore 💗

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved ® KALSTEIN France S. A. S., 2 Rue Jean Lantier • 75001 Paris • +33 1 78 95 87 89 / +33 6 80 76 07 10 • https://kalstein.eu KALSTEIN FRANCE, S. A. S



All rights reserved © KALSTEIN France S. A. S.,

www.kalstein.eu

KALSTEIN® Ultrasonic Cleaner Ultrasonic Washing Machine Portable Model YR05207



Product Description:

The Ultrasonic Washing Machine Portable is a microprocessor-controlled device that uses ultrasonic or high frequency sound waves to clean jewelry and other delicate items. Ultrasonic cleaning is based on the cavitation effect caused by high frequency ultrasonic wave vibration signals in fluid. Microscopic bubbles that expand and rapidly collapse cause cavitation, which creates an effective cleaning action on the surface of object being cleaned. Additionally, the bubbles are small enough to penetrate even microscopic crevices, cleaning them thoroughly and consistently.

Application:

- ✓ Professional Use: Medical and Dental Clinics, Tattoo Shops, Scientific Labs.
- Personal or Home Use: Jewelry: Earrings, Necklace, Rings, Bracelets and Diamonds.
- ✓ Glasses and Timepieces: Glasses, Sunglasses, Optical Lenses, Contact Lens Accessories, Watch Chains and Waterproof Watches.
- ✓ Commodities: Tattoo Guns and Tubes, Electric Shaver Heads, Razor Blades, Dentures, Combs and Toothbrushes.
- ✓ Stationery: Pen-heads, Printer-heads, Inkjet Cartridges and Seals.
- ✓ Metal Articles: Ancient Coins, Badges, Valves, Machine Nozzles, Electronics Components and Mechanical Parts.
- ✓ Metal Dishware: Forks, Knives, Spoons and Other Small Silverwares etc..



			 AC110-120V 60Hz
Name	Ultrasonic cleaner	Voltage	
			 AC220-240V 50Hz
Model	YR05207	Power	70W
Packing size	242*227*265mm	Ultrasonic frequency	40KHz
Internal groove size	150*135*100mm	Time setting	10-600S
Tank capacity	2L	Net weight	2.6KG



#Letsgivemore Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved ® KALSTEIN France S. A. S. 2 Rue Jean Lantier • 75001 Paris • +33 1 78 95 87 89 / +33 6 80 76 07 10 • https://kalstein.eu KALSTEIN FRANCE, S. A. S



www.kalstein.eu

All rights reserved © KALSTEIN France S. A. S.,





#LETSGIVEMORE

WITH THE ADCQUISITION OF A KALSTEIN EQUIPMENT

YOU MAKE A CONTRIBUTION TO:





https://onetreeplanted.org





ANY DOUBT? Contact us!



Articles Kalstein (Social Networks and Official Websites)





PARIS - FRANCE **HEADQUARTER** 2 Rue Jean Lantier 75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89 **Cel:** +33 (0) 6 80 76 07 10 sales@kalstein.eu


https://kalstein.eu/ https://www.kalstein.fr/

Kalstein Training Videos (Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...

23 views • 3 years ago

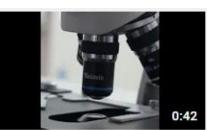

Kalstein's Microscopes 445 views • 1 week ago

