




LABORATORY LINES

REFRIGERATORS AND
FREEZERS



Kalstein®

A different accompaniment, at your service

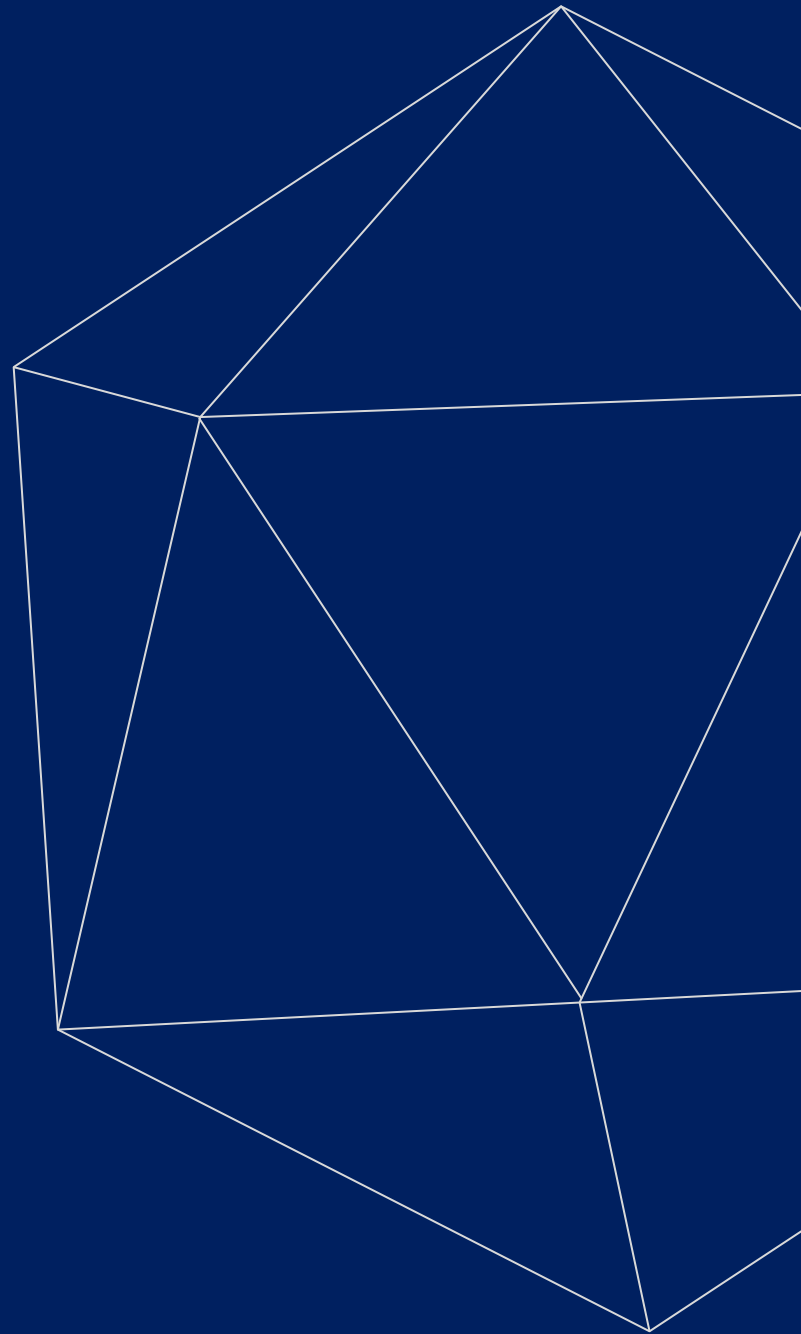



Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore ♥

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide



With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



-10 °C To -25°C Upright Freezer 528L

Serie YR05098

 Energy-saving and environmental friendly

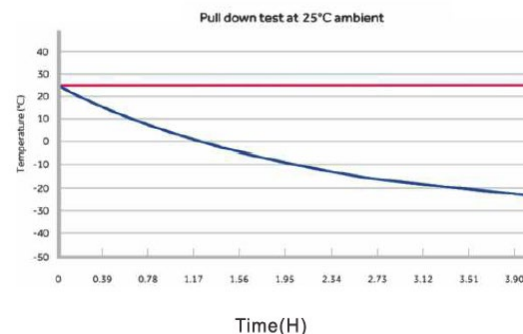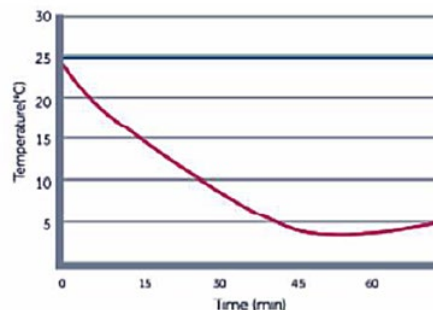

RELIABILITY

- Chemically stable, CFC-free, commercially available and environmentally safe refrigerant.
- High density CFC-Free urethane foam insulation for rigidity and stable storage temperature.
- Microprocessor control, digital display adjustable temperature range: -10C--25C.
- LED digital display for clear observation.



NEW RELEASE

- Freezer Room: -10°C-25°C /-30°C-40°C
- 2 Independent temperature display, 2 Independent door: heated foam glass door + foamed door
- 3 shelves in refrigerator room; 6 drawers in freezer room.



SAFETY

- Multiple malfunction alarms including high and low temperature, sensor error, power off, temperature controller malfunction, door ajar.
- Two types of alarm indications: audible buzzing and visible flashing light.
- Power failure protect: the device would be attack by the high current when all devices restart simultaneously, the start-delay function would protect device from the damage
- Controller password protect function to avoid change running parameter randomly



KEY FEATURES

- Independent inner door to avoid more cooling leaking out, and provide more stable inside temperature performance
- Upright side door with build-in handle for easy-opening; limit door hinge to prevent collision when door open
- Lockable caster design for customer to move as needed
- Wide voltage tolerance design with applicable voltage range of 187~242 V/AC.



#Letsgivingmore

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tit: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



-10 °C To -25°C Upright Freezer 528L

Serie YR05098



Microprocessor system, LED digital display, precise electronic temperature control and display, with accuracy of 0.1 °C; the temperature of refrigerator room and freezer room are displayed independently, which is convenient for observation. And the operation status and alarm type of refrigerator are displayed simultaneously.



2 independent door up and down: upper door is electric heating foamed glass door, below is foamed door.

The temperature of the refrigerator room is in the range of 2-8 °C, and freezer room is in the range of -10 °C ~ -25 °C \ -25°C ~ -40°C.

Two test holes are equipped on the right of the cabinet, installing the temperature recorder to detect the temperature change of each room

Independent refrigeration system is used in the refrigeration room and freezer respectively, with independent display and temperature control



#Let'sgive more ♥

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved (© KALSTEIN France S. A. S., Optimum Business Center 450 Rue Baden Powell, • 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710 / +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S




-10 °C To -25°C Upright Freezer 528L

Serie YR05098



Alarm	Alarm Triggering Condition
High Temperature	Temperature reaches the warm alarm limit
Low Temperature	Temperature reaches the low alarm limit
Power Failure	Equipment loses power
Door Ajar	Door open audible and visual alarm
Probe Failure	1. Main cabinet temperature control sensor fails 2. Sensor fails
Low Battery	Battery capacity runs low or battery switch is not turned on

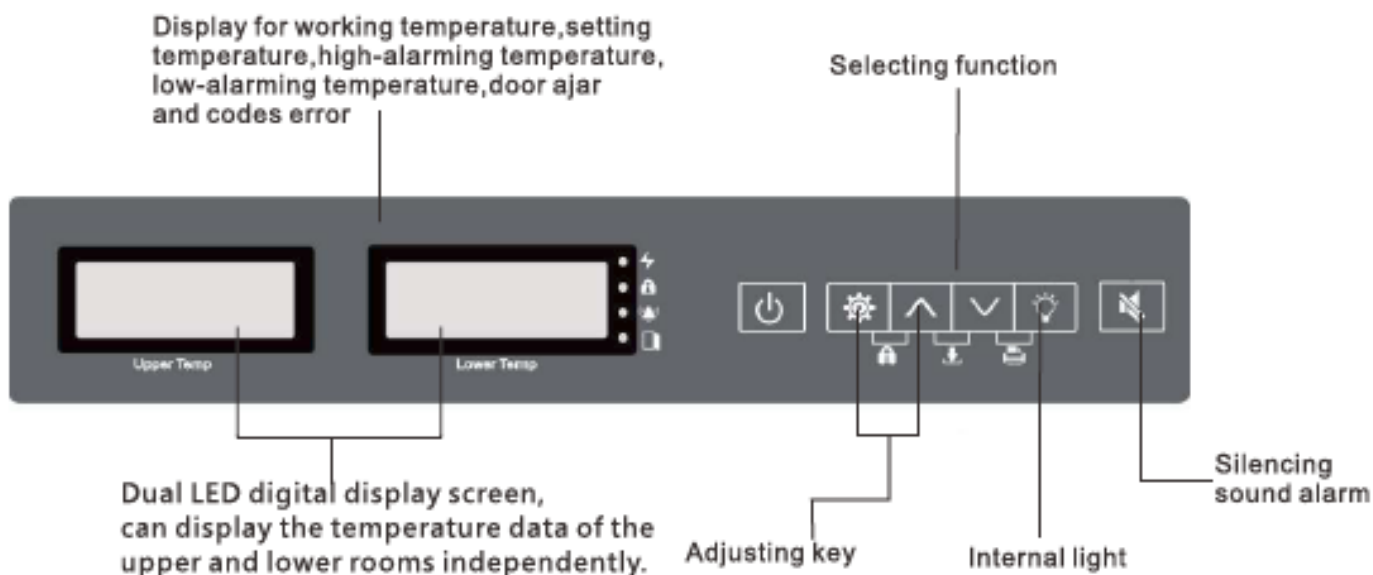

#Let'sgive more ♥

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



-10 °C To -25°C Upright Freezer 528L

Serie YR05098



Performance		Refrigeration System	
Cooling method	Direct cooling	Compressor brand/QTY	CHINA/2
Defrost	Manual	Condenser	Copper tube
Refrigerant	R290	Temperature controller	Microprocessor
Temperature range	-10~-25°C	Sensor	NTC
Ambient temperature	10~32°C	Display	Digital display
Electrical Data		Dimensions	
Voltage frequency(V/Hz)	220/50Hz	Capacity(L/cu.ft)	528/18.64
Power(W)	260W	Exterior size(W*D*H)(mm)	910*845*1875
Current(A)	1.9 A	Interior size(W*D*H)(mm)	(651*631.5*635)*2
Alarms		Other	
High/low temperature	Y	Inside door	9 Drawer×2
Power failure	Y	Test hole QTY/diameter	2/25MM
Sensor error	Y	Remote alarm port	Y
Temperature controller failure	Y	Castor	Y



#Letsgivingore ♥

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tit: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



Articles Kalstein
(Social Networks and Official Websites)



Kalstein Blog
Not just be informed, be fully informed

The 3 articles of the week!

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



1
Read this ARTICLE NOW

How does the climate affect your microscope?

microscope is an instrument that allows observing and measuring objects that are too small to be seen with the naked eye. This is achieved through a complex system of lenses and mirrors, which form

Kalstein Blog

Article of the week:
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW

What is the importance of the laboratory refrigerator?

refrigerator is a device that is one of the most common laboratory equipment. Its purpose is to maintain, in a controlled environment (refrigerated space), various fluids and substances

Kalstein Blog

Article of the week:
What is the importance of the laboratory refrigerator?


Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...

23 views • 3 years ago

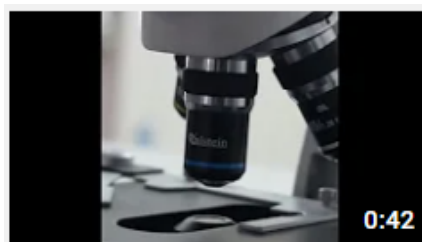

Kalstein's Microscopes

445 views • 1 week ago



MONTPELLIER - FRANCE
HEADQUARTER

Optimum Business Center,
450 Rue Baden Powell,
34000 Montpellier - France.

Fax: +33 (0) 4 67 15 88 49

Tif: +33 (0) 6 80 76 07 10

sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

Official Websites for Latin America.

